

NitraGrip XP Nitrile Exam Gloves

12" Powder-Free Non-Sterile Nitrile Examination Gloves

NitraGrip XP offer dependable strength you can count on and have a long 12" extended cuff that helps provide extra protection for your wrists and arms. Full textured for enhanced grip. Barrier tested against fentanyl exposure to help protect you when treating overdose patients.

Characteristics and Overview:

- When a high-risk situation requires the highest standard in protection.
- These elongated synthetic gloves provides extra protection and typically used for tasks involving prolonged exposure to blood, body fluids, surgical instrument cleaning solutions and other chemicals.
- Meets and Exceeds NIOSH recommendations for Fentanyl.
- Meets both ASTM D6978-05 and USP 800 Standards for use of chemotherapy gloves when handling hazardous drugs*




Item No	Description	Size	Packaging
MG50051	NitraGrip XP Nitrile PF Exam Glove	Small	100 eaches/box 10 boxes/case
MG50052	NitraGrip XP Nitrile PF Exam Glove	Medium	100 eaches/box 10 boxes/case
MG50053	NitraGrip XP Nitrile PF Exam Glove	Large	100 eaches/box 10 boxes/case
MG50054	NitraGrip XP Nitrile PF Exam Glove	X-Large	100 eaches/box 10 boxes/case
MG50055	NitraGrip XP Nitrile PF Exam Glove	2X-Large	90 eaches/box 10 boxes/case
MG50056	NitraGrip XP Nitrile PF Exam Glove	3X-Large	90 eaches/box 10 boxes/case

Glove Specifications:

Brand	NitraGrip XP
Material	Nitrile
Color	Blue
Sterility	Non-Sterile
Cuff	Beaded Cuff
Sizes	SM, MD, LG, XL, 2XL & 3XL
Accelerators	Yes
Chemo Tested	Yes
Glove Interior	Powder-Free
Glove Exterior	Full Textured
Length (inches)	12 in
Palm Thickness (Avg)	5.0 mil
Finger Thickness (Avg)	6.2 mil
Tensile Strength (Mpa)	38 / 35 Mpa
Ultimate Elongation (%)	650 / 600
Pin Hole (AQL)	1.5

Quality Standards:

Manufactured in accordance with Quality System ISO 9001.
Exceeds and meets current ASTM Standards:
- D6319, D5151, D1671 & D6978.

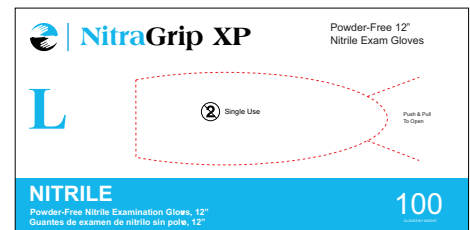

*Summary of Chemotherapy Testing Reports - ASTM D6978-05

Drug Name	Average Breakthrough Detection Time
Carboplatin	No breakthrough up to 240 min
Cisplatin	No breakthrough up to 240 min
Cyclophosphamide	No breakthrough up to 240 min
Dacarbazine	No breakthrough up to 240 min
Doxorubicin	No breakthrough up to 240 min
Etoposide	No breakthrough up to 240 min
Fluorouracil	No breakthrough up to 240 min
Ifosfamide	No breakthrough up to 240 min
Methotrexate	No breakthrough up to 240 min
Mitomycin C	No breakthrough up to 240 min
Mitoxantrone	No breakthrough up to 240 min
Paclitaxel	No breakthrough up to 240 min
Vincristine Sulfate	No breakthrough up to 240 min
Carmustine	Not for use with Carmustine
ThioTepa	Not for use with ThioTepa

